Mirror Chart MC-5/MC-5S



MC-5/MC-5S

MC-5/MC-5S supports you in the path to adapting quickly to new social distance recommendations whilst providing sophisticated visual acuity (VA) and space saving.



Installation distance: 40cm*1



Installation depth: 21cm



MC-5S

NEW refraction system with "Circular polarization"

Circular polarization

The built-in circular polarization panel provides stability even if an individual moves their head. This is due to the stability provided by the binocular balance visual function test with a polarization chart.*2*3





OK!

3D contents can also be performed on the MC-5S.

MC-5




NEW LCD Remote controller!

LR-5

The new LCD remote control is more comfortable to hold due its new built-in LR-5 alignment emmiter.




- *1 Minimum installation distance 40cm with MC-5, 60cm with MC-5S.
- *2 With wearing the circulor polarization glasses while the test.
- *3 The Circular polarization glasses are not a standard accessory while wearing circular polarization glasses.
- *4 Compared with LR-4.

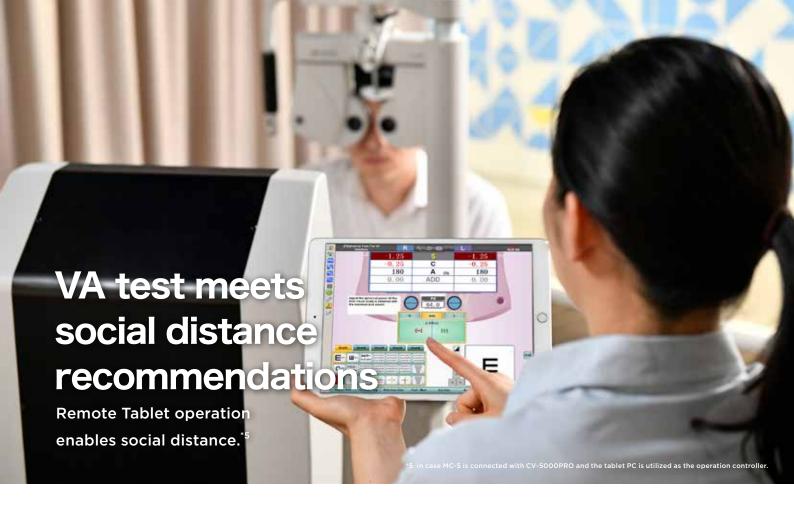

Space saving VA test and flexible installation

- The VA test with a small installation distance.
- Flexible placement positioning in your practice or clinic.
- · The multiple refraction lanes in a small area.

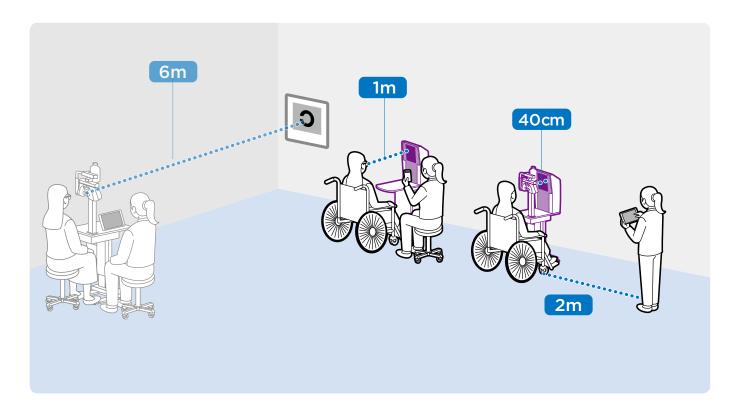

Optional accessory for MC-5



LCD remote controller

LR-5

Combination devices with MC-5S



CV-5000PRO



1 dial controller

KB-50S

Specifications

Installation distance	MC-5: 40cm, 65cm, 85cm, 100cm, 105cm, 155cm MC-5S: 60cm, 80cm, 95cm, 110cm, 150cm* *Only for H1 type					
VA value	$\begin{array}{ll} \text{MC-5}: & 0.03\sim2.0^{**}\\ \text{MC-5S}: & 0.03\sim2.0 \end{array}$ **0.03 can be displayed when MC-5 is connected with CV-5000PRO.					
Test chart	MC-5: Landolt ring, numerals tumbling E, R&G, astigmatism,.etc. MC-5S: Landolt ring, numerals, alphabets, tumbling E, R&G, astigmatism,.etc.					
Dimensions	450mm (W), 207mm (D), 621mm (H)					
Weight	MC-5: 13.6kg, MC-5S: 13.1kg					
Power supply	AC: 100V-240V, 50-60Hz					

Chart sample

(F type for MC-5)

0.05	000 Q1	COOOO0.3	OOOO0.6	000000.8 000000.9 000001.0	COOCO 1.5	4 5 8 0.1 2 9 3 0.2	543670.3 682940.4 328650.5	934270.7	94527 _{1.0} 68349 _{1.2} 59283 _{1.5}
∃Ш П 0.1 ПЕШ 0.16	ЕМЭШЕ 0,2 ШЕМЭМ 0,3 ЕЭШМЕ 0,4	3EWMW 0.5 W3ME3 0.6 MWE3M 0.7	3MEWE 0.8 MW3EM 0.9 WEM3E 1.0	ЭШЕМШ 1.2 ЕЭШЕМ 1.5 ШЕМШЭ 2.0	•	5 © Ø 5 34 43 290 092	1	\$ \$\pi \pi \rightarrow \Vartial 0.4 \$\rightarrow \Vartial \pi \rightarrow \Vartial \pi \rightarrow 0.5 \$\pi \Vartial \rightarrow \Vartial \Vartial 0.6	
	1112 1 10 1 2 9 3 8 7 6 4 7 6 5	00000	00000	0000	ii Ii			(+++

Manufacturer

III TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, JAPAN. Phone: +81-(0)3-3558-2522/2502 Fax: +81-(0)3-3965-6898 www.topcon.co.jp

TOPCON SINGAPORE MEDICAL PTE. LTD.

1 JALAN KILANG TIMOR #09-01 PACIFIC TECH CENTRE SINGAPORE 159303 Phone:+65-68720606 Fax:+65-67736150 E-mail:medical_sales@topcon.com.sg www.topcon.com.sg

TOPCON INSTRUMENTS (MALAYSIA) SDN. BHD.

No.6, Jalan Pensyarah U1/28, Hicom Glenmarie Industrial Park 40150 Shah Alam, Selangor, MALAYSIA Phone: +60-(0)3-50223688 Fax: +60-(0)3-50313968

TOPCON INSTRUMENTS (THAILAND) CO., LTD.

77/162 Sinnsathorn Tower, 37th Floor, Krungthonburi Rd., Klongtonsai, Klongsarn, Bangkok 10600, THAILAND Phone: +66(0)2-440-1152-7 Fax: +66-(0)2-440-1158

MEHRA EYETECH PTE. LTD.

801 B Wing, Lotus Corporate Park, Graham Firth Steel Compound Goregaon (East) Mumbai 400063 Maharashtra, INDIA Phone: +91-22-61285455 www.mehraeyetech.in

TOPCON (BEIJING) MEDICAL TECHNOLOGY CO., LTD.

Room 2808, Tower C, JinChangAn Building, No.82, Middle Section of East 4th Ring Road, Chaoyang District, Beijing 100124, P.R. CHINA Phone: +86-10-8794-5176

IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

Not available for sale in all countries. Please check with your local distributor for availability in your country. MC-5S is an accessory of CV-5000PRO.

Compu Vision CV-5000PRO ("PRO" is the brochure name.)





